



WE'RE ALL ABOUT PEOPLE.

Dornier *Compact Sigma*

Modular lithotripter



Dornier *MedTech*

Dornier Compact

The Perfect Combination:
The Modular Dornier *Compact Sigma*



+++ Extracorporeal Shock Wave Lithotripsy +++ Orthopedics +++

The Dornier *Compact Sigma* The Complete Solution

Dornier *Relax+* Patient Table

The fully motorized, four-axis table is ideal for both Extracorporeal Shock Wave Lithotripsy (ESWL) and endourology. It incorporates an exchangeable radiolucent stretcher that offers a wide range of movement as well as Trendelenburg tilting capability. Extensive urological accessory options are also available.

TLX Patient Table

The cost-effective motorized TLX table has been optimized for ESWL applications with the Dornier *Compact Sigma*. It features smooth operation and comfortable patient positioning.

Ultrasound Scanners from BK Medical

Select one of the advanced scanners from the ultrasound leader in urology, BK Medical. Both the Merlin and Pro Focus provide special transducers for shockwave applications, prostate scanning, and abdominal diagnosis.

Surgical C-arms

Selected surgical C-arms, from a variety of manufacturers, can be mechanically linked to the *Compact Sigma* to create a single stable unit. This provides high quality imaging for lithotripsy as well as other applications.



Endourology +++ Sonography +++ Fluoroscopy +++ Brachytherapy +++ TUR +++ URS

Depth of Experience



Proven Therapeutic Results with the Dornier *EMSE 140f*

The Dornier *EMSE 140f* is the result of vast experience and expertise in shock wave technology. Real world therapeutic effectiveness and minimization of side effects were the design objectives. Today hundreds of customers worldwide trust the Dornier technology. With high penetration depth and a wide dynamic range, the *EMSE 140f* is ideal for a spectrum of applications – from low-energy analgesic therapy to low-anesthesia stone treatment to high-energy pseudarthrosis application.



Compact Sigma

Designed Around the Focus



Isocentric – Designed for Precision and Flexibility

The Dornier *Compact Sigma*'s isocentric design allows both the shock wave and the imaging systems to revolve around a single focal point. In practical terms, this translates into the precise alignment of the targeting system with the shock wave focus. The projection angles of the X-ray localization and the ultrasound windows can be varied over a wide range for optimal imaging, without losing the relationship to the therapy focus. Thanks to this isocentric design, elaborate and costly navigation systems for stone localization are unnecessary with the *Compact Sigma*.

Flexible Shock Wave – Promotes Patient Comfort


Depending on the indication, shock waves can be introduced from different orientations. This flexibility allows all treatments to be conducted with patients positioned comfortably on their backs, a unique feature of Dornier's *Compact* line of lithotripters.



Isocentric

Dornier Compact

Ultimate Connection

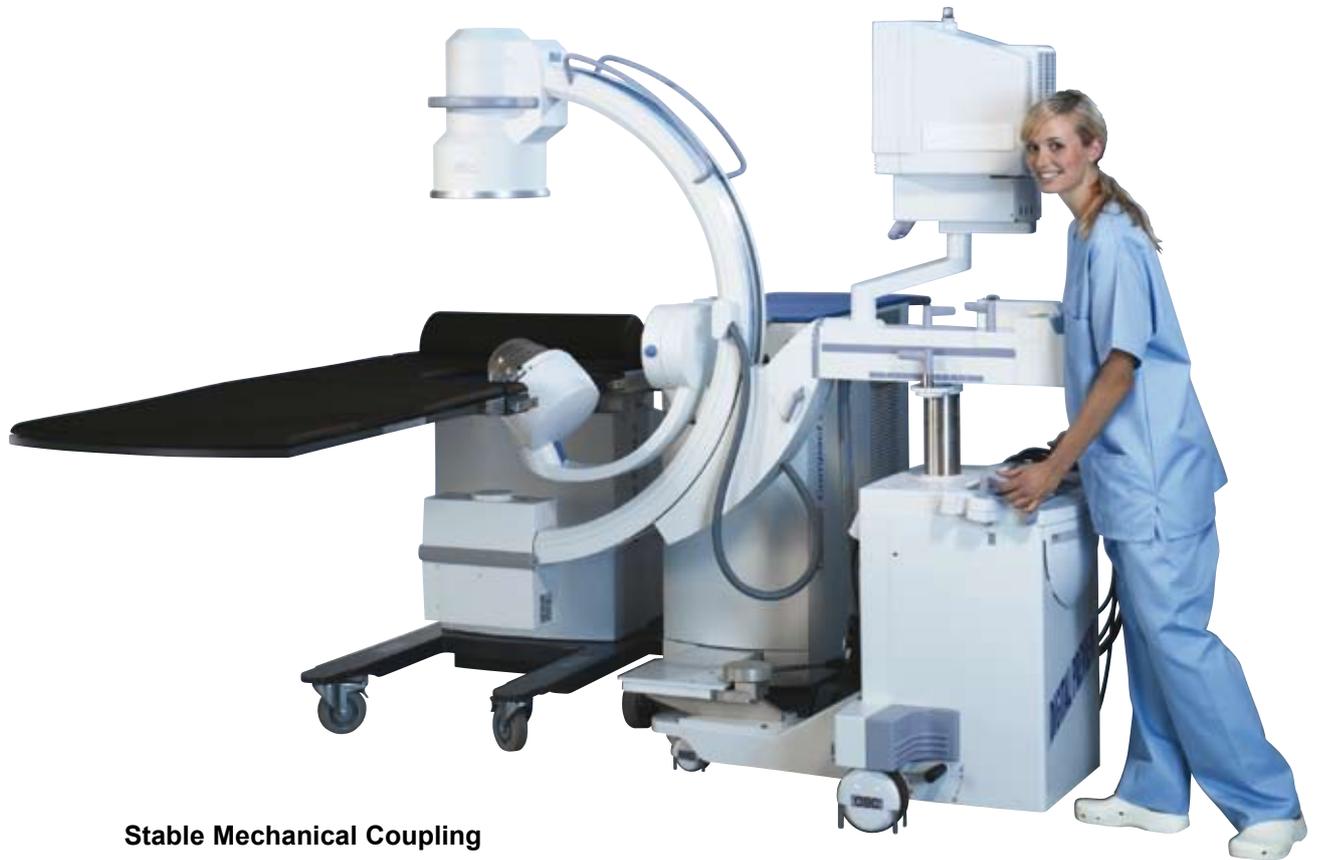

Stable Mechanical Coupling

Simply slide the X-ray C-arm under the connecting plate and turn the lever. That's it. The mechanical coupling is fast to set up, easy to operate and extremely stable. A variety of C-arms from well-known manufacturers, are supported.



Sigma



Laser-guided Alignment

QuickLinx®, Dornier's optical position control, ensures proper alignment between the C-arm and *Compact Sigma*. The laser aiming device is mounted to the top of the therapy arm and is equipped with three manually adjustable beams, which can be turned on or off with the touch of a button.

**Quick
LinX**

Precise Alignment

Dornier Compact

Patented Localization – Flexible and Precise

Ultrasound Freedom

The use of the isocentric localization arm allows the shock wave and ultrasound to remain independent of each other. This is advantageous in situations where interaction between the shock waves and ultrasound picture can create issues. The results are improvements in the effectiveness of the shock wave as well as the lifetime of the electronic probe. During therapy, the ultrasound transducer is always in direct contact with the skin – to provide optimal image quality.




Stone Precision

The localization arm provides the flexibility needed for excellent image quality and accurate targeting.



Ortho Expertise

The *Compact Sigma* can also be used to treat a wide variety of tendonopathies. The flexible therapy head supports simple targeting for various indications.

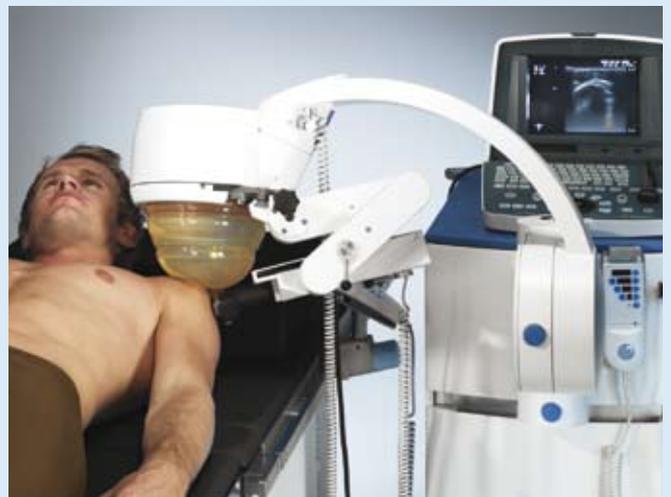

Sigma



FarSight

FarSight Revolutionary "See What You Shoot Targeting"

The *Compact Sigma* offers an additional enhancement for ultrasound imaging flexibility with FarSight. This state-of-the-art integrated ultrasound technology ensures precise stone targeting and real-time monitoring of the shock wave path and the disintegration process. FarSight revolutionizes conventional inline ultrasound by delivering TRUE "see what you shoot" imaging while you shoot along the shock wave path.



FarSight integrates a specially designed curved array transducer into the therapy head of the *Compact Sigma*. The integrated transducer allows for deep imaging while remaining outside the shockwave path at all times. The motorized transducer can be adjusted at the touch of a button to rotate along the longitudinal scanning plane.



Picture Perfect

Dornier Compact

Ideal Orientation



Optimized Indication Settings for Shock Wave and Localization

Proximal ureteral stones can be localized by both vertical and oblique X-ray projections providing posterior lateral coupling for treatment. This leaves the X-ray image clear and unobstructed. Avoidance of bone or intestinal gas is simple thanks to the flexible imaging design.



Lithotripsy in the medial and distal ureter is typically performed ventrally due to the osseous pelvis. The overtable position of the therapeutic unit provides excellent access while patients remain comfortably on their backs.




Dual Imaging

The simultaneous use of X-ray and ultrasound can be used for stones that are not clearly seen with ultrasound alone. X-ray localization allows instant positioning while ultrasound provides real time information about respiration and disintegration.

Sigma

Tri-mode Imaging

Dornier's unique tri-mode imaging extends your imaging capabilities far beyond what conventional lithotripters can offer. Tri-mode imaging not only provides you with the most appropriate imaging modality for any stone location but also allows for optimum monitoring during shock wave application.

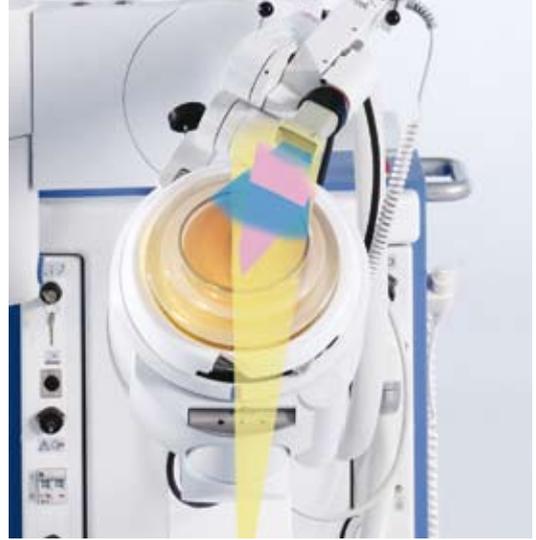

Imaging Versatility

UIMS Urology Information Management System




Dornier UIMS Viewing Station SC

UIMS – Dornier’s comprehensive Urology Information Management System is optionally available on the **Compact Sigma** for viewing, processing, and archiving of images and for patient data management. Furthermore, it offers DICOM capabilities and can be connected to a PACS network to maximize patient data collection and workflow.


Viewing Station SC

The **UIMS** Viewing Station can operate as a stand alone viewer, when attached to a hospital network or as a secondary capture system for a variety of imaging devices including C-arms, ultrasound systems and endoscopes. Based upon the viewing station model selected, the secondary capture feature can directly acquire live images and store single images or sequences from a standard line or highline image device. The acquired images can be transmitted to the facility’s PACS system using the DICOM interface or viewed on the LCD monitors.

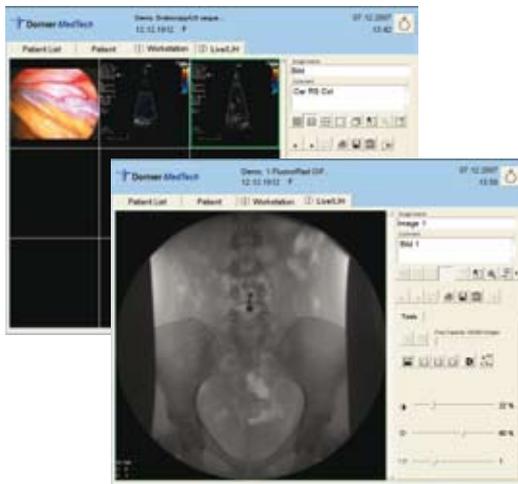


System Flexibility

UIMS is configured to accept a wide range of imaging devices. For easy set-up, a wide variety of imaging device templates is included in the system.

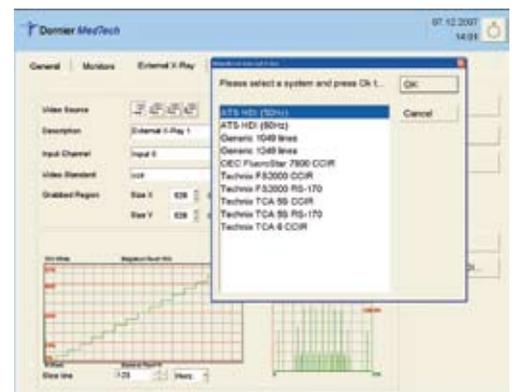


Compact Sigma

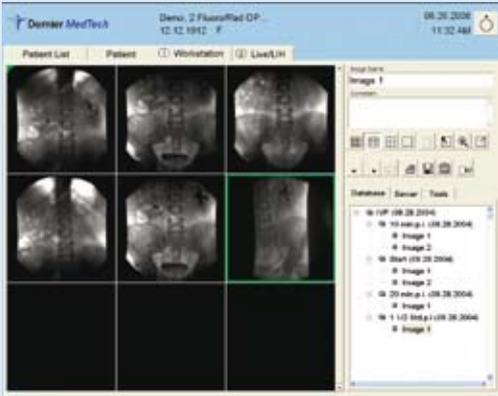

Increased Workflow

UIMS features a user-friendly, intuitive graphical interface allowing patient data and treatment information to be entered quickly. Using a simple search function or DICOM Worklist, patient information can be retrieved easily enhancing workflow. Additionally, **UIMS** is available on Dornier's urology imaging table and most Dornier lithotripters providing inter-operability on Dornier **UIMS** systems.

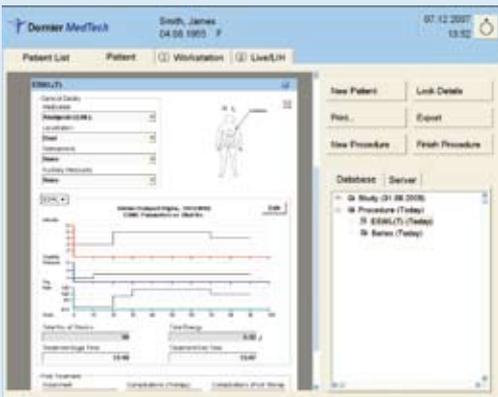

Advanced Patient Data Management

ESWL treatment information can be documented in the Advanced Patient Data Management System for a complete patient treatment history. Pre-treatment, ESWL data and post-treatment information can be included and select data presented in a graph or histogram. User-defined treatments reports can be created to meet a variety of needs.



Image Archive and Print

Images can be stored to an assortment of media storage devices including PACS, DVD, CD, memory stick in a variety of formats. When exporting the images to a local or network printer, the operator can use the Print Preview function to select the desired images.

Ergonomic by Design

Quick & Easy

Settings are fast and reliable. In a matter of seconds the *Compact Sigma* can be adjusted for various therapeutic positions.



On the Move

Large and stable wheels provide easy transportation. Door thresholds, uneven pavement and slight inclines don't cause transportation challenges for the *Compact Sigma*.

Safely Parked

The foot-activated safety brake ensures the *Compact Sigma* stays parked.



Compact Sigma



Tucked Away

The connecting plate for the mechanical linkage of the C-arm, tucks away neatly. Just fold it out when needed.

Cable Free

The *Compact Sigma's* cables are easily recoiled and hidden in a specially designed side-compartment.



Additional Storage

The *Compact Sigma's* sleek ergonomic design includes a convenient pull-out storage drawer for a variety of accessories.

Everything in Place

Headquarters

Dornier MedTech GmbH
Argelsrieder Feld 7
D-82234 Wessling
Germany

America

Dornier MedTech America, Inc.
1155 Roberts Boulevard
Suite 100
Kennesaw, GA 30144
USA
Phone: +1-770 426-1315
Fax: +1-770 426 6115
Email: info@dornier.com

Asia/Pacific

Dornier MedTech Asia Pte. Ltd.
2 Corporation Road
#05-11/12 Corporation Place
Singapore 618494
Phone: +65-6665-5788
Fax: +65-6665 6711
Email: infoasia@dornier.com

Japan

Dornier MedTech Japan
Meguro Estate Building
3-8-5 Kami-Ohsaki
Shinagawa-ku, Tokyo 141-0021
Japan
Phone: +81-3-3280 3550
Fax: +81-3-3280 3555
Email: info@domedtech.co.jp

Europe

Dornier MedTech Europe GmbH
Argelsrieder Feld 7
D-82234 Wessling
Germany
Phone: +49-8153-888 625
Fax: +49-8153-888 444
Email: infoeurope@dornier.com

Dornier MedTech Italia s.r.l.
Via Arrigo Cavaglieri N. 26
00173 Rome
Italy
Phone: +39-06-72 35 289
Fax: +39-06-72 35 302

Dornier MedTech Espana
Pedro Y Pons 9-11
Edificio Masters
E-08034 Barcelona
Spain
Phone: +34-93-203 9316
Fax: +34-93-205 3814

Dornier MedTech France SARL
Ferme de St Paul
1545 route d'Epagny
F-74330 Sillingy
France
Phone: +33-450-22 18 94
Fax: +33-450-24 24 63

Dornier MedTech Europe GmbH
Moscow Representative Office
Uliza Mytnaya 3
Office 21
119049 Moskau
Russia
Phone: +7-495-230 29 94
Fax: +7-495-739 51 33



CE 1275